


Cochlear™ Nucleus® Profile™ Plus with Slim Modiolar Electrode (CI632)

MRI AT 1.5 T AND 3.0 T WITH MAGNET IN PLACE

- 1 The thinnest implant body with no pedestal¹, designed to minimise bone excavation and skin protrusion.
- 2 Implant coil enabling telemetry.
- 3 Smooth external geometry to minimise biofilm formation and reduce risk of infection.²
- 4 Titanium casing, for impact resistance.
- 5 Symmetrical side-by-side exit leads from main casing, for easier surgery. Same procedure for left and right ear.
- 6 Removable magnet for MRI safety to minimise image distortion. MRI at 1.5 Tesla and 3.0 Tesla with magnet in place.³

Circle on magnet indicates the side that should be away from the bone.


Sterilised replacement magnet cassette (P782485) and non-magnetic cassette (P782484) to assist MRI procedures available from Cochlear.

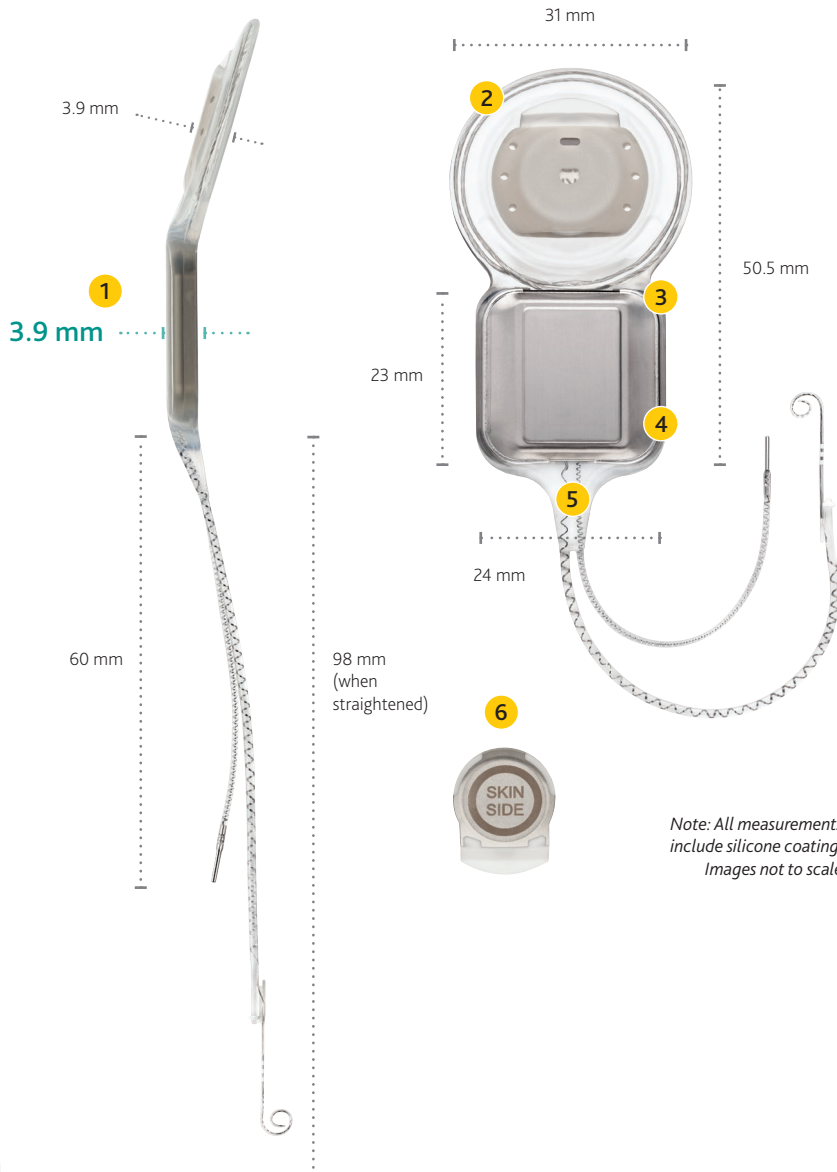

Weight (without sheath)	9.2 g including electrode array.
Impact	Resistant against external impact up to 2.5 joules. ⁴

SLIM MODIOLAR ELECTRODE

- 1 Intracochlear electrode with 22 platinum electrode contacts spread over 14 mm active length.
- 2 Distance from electrode tip to most proximal electrode contact - 14.4 mm.
- 3 Dimensions at basal end 0.475 x 0.5 mm.
- 4 Dimensions at apical end 0.35 x 0.4 mm.
- 5 Three white insertion depth markers, visible only after sheath is removed. Distance between white markers - 1 mm.
- 6 Distance from electrode tip to most proximal white marker - 18.4 mm.
- 7 White alignment marker on electrode - diameter 2.1 mm.
- 8 Position where the sheath can be reloaded onto the electrode.
- 9 Proximal sheath end - diameter 0.77 mm.
- 10 Internal sheath - length 5.5 mm.
- 11 Sheath stopper - diameter 1.5 mm.
- 12 White alignment marker on sheath - diameter 1.45 mm. When electrode is fully inserted, the markers on the electrode and the sheath are aligned.
- 13 Sheath guide tube (including internal sheath length) - length 28.9 mm.
- 14 Sheath handle - length 9.5 mm.
- 15 Cochleostomy Sizing Tool - length from tip to stopper 2 mm.
- 16 Cochleostomy Sizing Tool tip - diameter 0.8 mm.
- 17 Cochleostomy Sizing Tool stopper - diameter 1.4 mm.

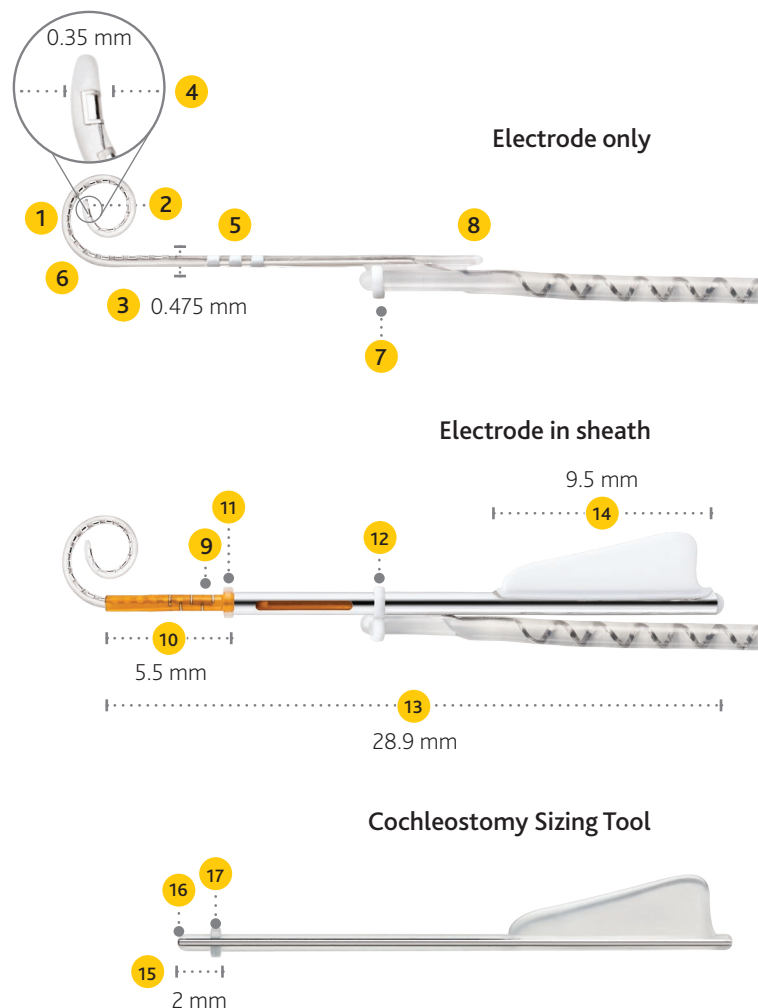

- ☒ The Slim Modiolar Electrode is suitable for round window or cochleostomy surgeries.
- ☒ A Sterile Silicone Implant Template and Cochleostomy Sizing Tool are provided in the surgical kit.

CUSTOM DESIGN MICROELECTRIC PLATFORM

Power efficient

- Stimulus amplitude range 0 μ A to 1750 μ A nominal at 37°C.
- Stimulation rates up to 31.5 kHz.

Implant ID and type check

- Enables the sound processor to confirm whether it is coupled to the nominated implant.

Stimulation modes

- Monopolar, bipolar or common ground stimulation modes using biphasic current pulses, designed for flexible programming options.

Telemetry capability

- Includes fully integrated electrophysiology telemetry modes - NRT®, AutoNRT® ESRT, ABR, CEP and intraoperative NRT.

1 Compared to all currently available receiver stimulators available from Cochlear and other cochlear implant manufacturers. Based on published device specification information.

2 James G A, Boegli L, Hancock J, Bowersock L, Parker A, Kinney B M, Bacterial Adhesion and Biofilm Formation on Textured Breast Implant Shell Materials, Aesth Plast Surg, October 2018; <https://doi.org/10.1007/s00266-018-1234-7>

3 MRI compatibility may vary by country depending on regulatory approvals in each country. Please check the MRI guidance provided in your country by contacting your local Cochlear representative or clinic before proceeding with an MRI scan.

4 EN 45502-2-3 Active implantable medical devices - Part 2-3: Particular requirements for cochlear and auditory brainstem implant systems.